

## HAND SPLINT FOR BURNS

WRIST IN NEUTRAL POSITION\*\* M.P. JOINTS APPROXIMATELY 80 FLEXION  
DESIGNED TO ACCOMMODATE BULKY DRESSINGS

### INSTRUCTIONS


Carefully place bandaged hand in splint . Evaluate fit. Edges should not “dig in”.

Kydex plastic can be modified with heat gun. Remove from patient. Carefully heat and reshape.

Allow to completely cool before using.

Use a soft bandage to secure splint to patient.

Clean plastic with warm water and mild soap. Rinse thoroughly. Towel and air dry completely before reapplying.

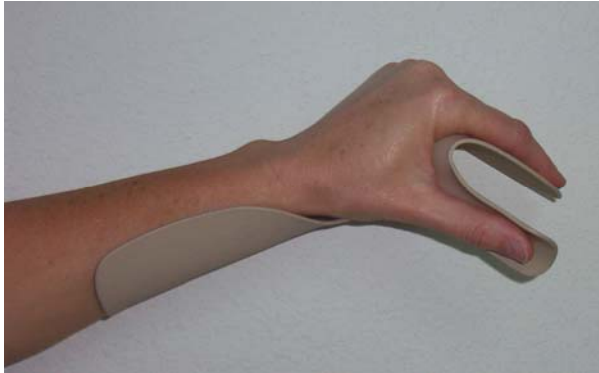

SIZE	CATALOG #	WIDTH AT M.P.
SMALL RIGHT	562847	2 7/8” - 3 1/8”
SMALL LEFT	562846	2 7/8” - 3 1/8”
MEDIUM RIGHT	562849	3 1/8” - 3 1/2”
MEDIUM LEFT	562848	3 1/8” - 3 1/2”
LARGE RIGHT	562851	3 1/2” - 3 7/8”
LARGE LEFT	562850	3 1/2” - 3 7/8”
X-LARGE RIGHT	562853	3 7/8” - 4 1/4”
X-LARGE LEFT	562852	3 7/8” - 4 1/4”

Single patient use. To be used only when ordered by a physician and applied by qualified medical professional



## HAND SPLINT FOR BURNS

WRIST IN NEUTRAL POSITION\*\* M.P. JOINTS APPROXIMATELY 80 FLEXION  
DESIGNED TO ACCOMMODATE BULKY DRESSINGS

### INSTRUCTIONS

Carefully place bandaged hand in splint . Evaluate fit. Edges should not “dig in”.

Kydex plastic can be modified with heat gun. Remove from patient. Carefully heat and reshape.

Allow to completely cool before using.

Use a soft bandage to secure splint to patient.

Clean plastic with warm water and mild soap. Rinse thoroughly. Towel and air dry completely before reapplying.



SIZE	CATALOG #	WIDTH AT M.P.
SMALL RIGHT	562847	2 7/8” - 3 1/8”
SMALL LEFT	562846	2 7/8” - 3 1/8”
MEDIUM RIGHT	562849	3 1/8” - 3 1/2”
MEDIUM LEFT	562848	3 1/8” - 3 1/2”
LARGE RIGHT	562851	3 1/2” - 3 7/8”
LARGE LEFT	562850	3 1/2” - 3 7/8”
X-LARGE RIGHT	562853	3 7/8” - 4 1/4”
X-LARGE LEFT	562852	3 7/8” - 4 1/4”

Single patient use. To be used only when ordered by a physician and applied by qualified medical professional

