


MedVue[®]

Medical Weight Analyzer

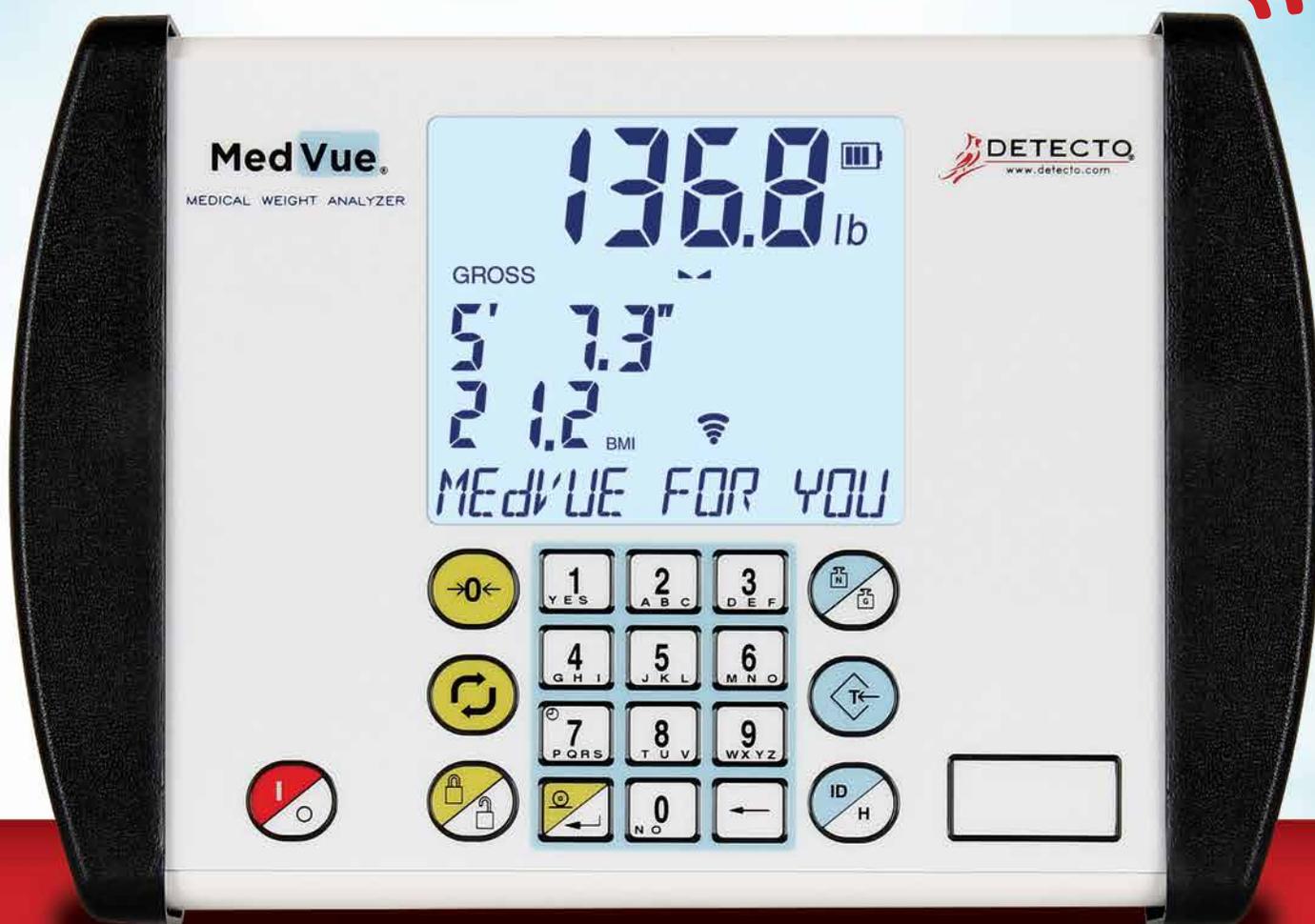

Detecto Model MV1 Medical Weight Analyzer

- Large Bold Display Readouts
- 350-Degree Swivel Mounting Bracket
- Wi-Fi Option for EMR/EHR
- 200 Stored Transactions
- 99 Stored Tare Weights
- 14-Digit Patient IDs
- Alphanumeric Keypad



Versatile Connectivity

DETECTO's MedVue medical weight analyzer offers a host of different options, both wired and optional wireless, for DETECTO scales to connect with multiple devices for efficient EMR/EHR.

Connect with:

- Printers
- PCs
- Networks
- Tablets
- Bar Code Scanners
- Even Mobile Phones



EMR/EHR connectivity requires custom coding by software package developers.

BENEFITS OF WI-FI CONNECTIVITY

DETECTO's MedVue offers optional Wi-Fi connectivity with up to 20 ft. range (varies depending upon environmental conditions). Wi-Fi offers a number of benefits over Blue-tooth, including greater range and security.



Digital Height Rod Port

Wireless (optional)

RS232 Serial

USB

Ethernet (optional)



P50 printer

Almost every MedVue port may be used simultaneously (with the exception of Ethernet and Wi-Fi connections). For instance, a printer could be connected to the serial port, PC connected to the USB port, while a tablet is connected via Wi-Fi all at the same time. DETECTO's model P50 printer offers plug-and-play compatibility with the MedVue.

350-Degree Swivel Mounting Bracket



350-degree rotation in all directions!

Enjoy freedom of movement wherever you need to position the display with the MedVue's 350-degree tilt-and-swivel mounting bracket.




Easy-to-read,
7/8-inch (22-mm)
high LCD weight digits



Weight, height, and BMI shown
onscreen simultaneously

Large, bold 3.5-inch-wide x
3-inch-high (89-mm-wide x
76-mm-high) backlit display

Wi-Fi Option

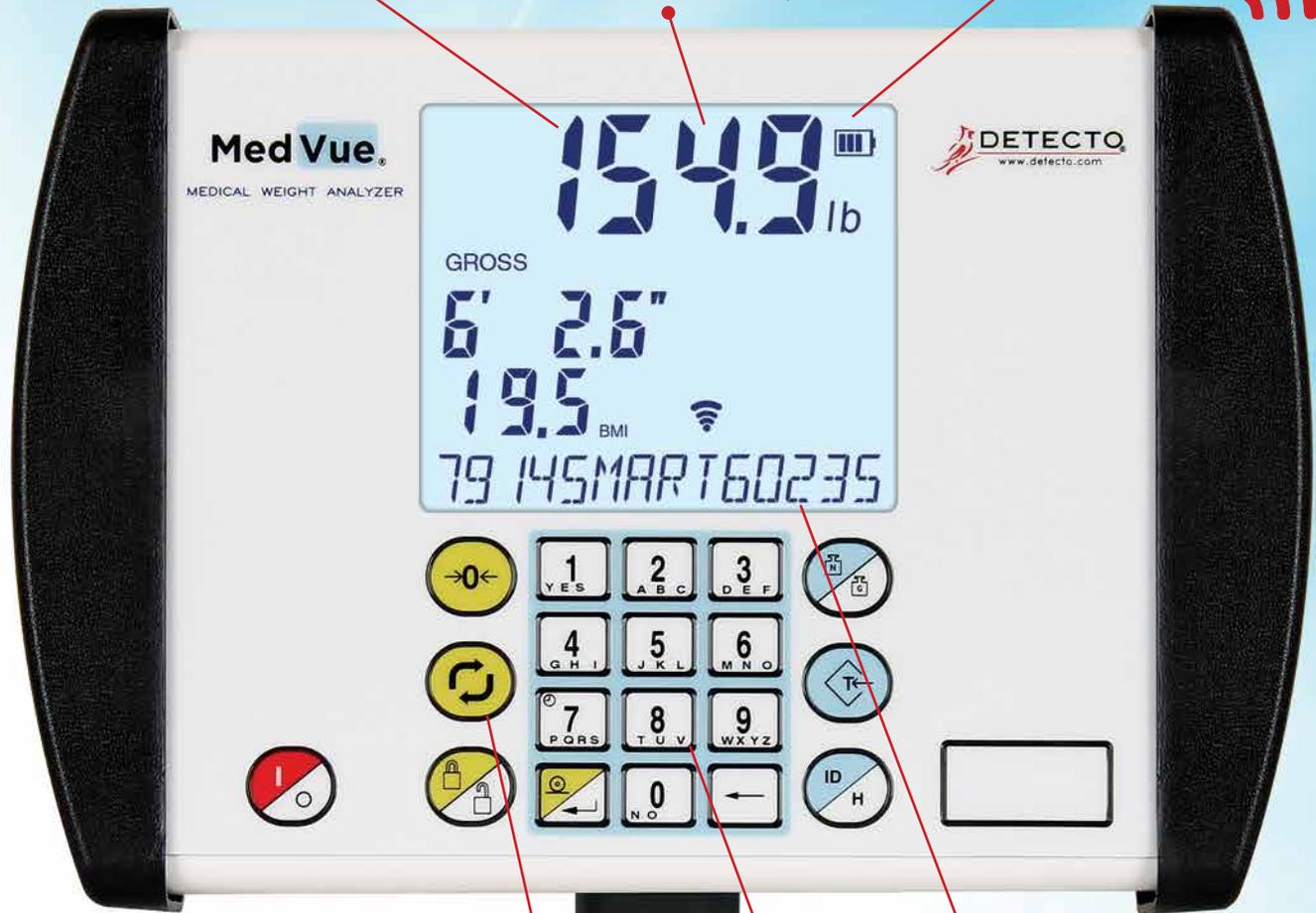

Digitally filters the weight signal to produce
a stable weight display despite movement

Convenient
units switching

Alphanumeric
keypad

14-digit
patient ID

Selectable sleep mode
conserves battery power

The MedVue is Standard on Select DETECTO Clinical Scales

EASY SET-UP NAVIGATION

DETECTO's MedVue indicator features a tiered menu system with up/down navigation controls to quickly update any settings you wish to change, such as converting to kilograms-only weighing or lowering the backlight intensity to make the batteries last longer.



MULTIPLE LANGUAGE SETTINGS

The MedVue features user-prompting in four different languages:

- English
- Spanish
- French
- German



Standard on select
eye-level, waist-high,
wheelchair, and
bariatric scales.



Model
6437DHR

Model
6857DHR

Model 6550

MedVue[®]

Medical Weight Analyzer

Detecto's MedVue operates on 6 C batteries (not included) or optional MV1PWR AC adapter.



SPECIFICATIONS

Power Requirements	6 "C" cell Alkaline, Ni-Cad or NiMH batteries (not included) OR an optional 12 VDC 1.25A wall plug-in AC power adapter (Detecto part number MV1PWR)
Display	Overall Display Size: 3.5 inch W x 3 inch H (89 mm W x 76 mm H) backlit blue LCD Weight: 5 digit, seven segment, 0.875 inch (22.23 mm) high LCD Height: 4 digit, seven segment, 0.5 inch (12.7 mm) high LCD BMI: 3 digit, seven segment, 0.5 inch (12.7 mm) high LCD ID: 14 digit, 14 segment, 0.35 inch (8.89 mm) high LCD
Battery Life	7,560 weighments (with backlight fully ON) / 10,620 weighments (with backlight fully OFF)
Sensitivity	0.5 uV/division (0 to 3.0 mV/V), Class III
Operating Temperature	14 to 104 °F (-10 to +40 °C)
Signal Input Range	0 to 15mV max.
Transducer Excitation	5.0 VDC
Number of Load Cells	up to 4 each 350 Ω
Load Cell Cable Length	30 feet (9 meters) max.
Capacities	Weight: 10,000 divisions commercial - Height: 0 - 9' 11.9" - BMI: up to 999
Resolution	10,000 divisions
NTEP	Certificate of Conformance Number 13-011
Division Value	1, 2 or 5 x 1, 0.1, 0.01, 0.001
Sample Rate	1 to 100 samples per second selectable
Auto Zero Range	0.5 or 1 through 9 divisions
Weighing Units	Pounds (lb), Pounds/Kilograms (lb/kg), Kilograms (kg) and Kilograms/Pounds (kg/lb)
Keyboard	Membrane type with 19 keys
Enclosure Size	9.8 inch W x 6.8 inch H x 1.8 inch D (248.9 mm W x 172.7 mm H x 45.7 mm D)
Construction	Painted steel



STANDARD FEATURES:

- Metric Conversion • Auto Shut-Off • Selectable Sleep-Mode • Selectable Filtering • 99 Tare Storage
- Serial I/O (Input/Output) • USB Device Port • Height Rod Port (for Optional Detecto Digital Height Rod) • 200 Transactions Stored

OPTIONAL FEATURES:

- MVIP Ethernet Option • MVWIFI Wireless Option • Detecto DHR (Digital Height Rod) • MV1PWR, 12VDC AC Power Adapter
- Welch Allyn[®] WACONNECT **WelchAllyn[®]**

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.

 **DETECTO[®]**
A Division of Cardinal Scale Manufacturing Co.
203 E. Daugherty, Webb City, MO 64870 USA
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001
www.Detecto.com

Mobile: m.detecto.com

SOLD BY:

